



MED MIZER

INNOVATIVE SOLUTIONS  
FOR HOME CARE

# MED-MIZER HOME CARE SOLUTIONS

We understand that keeping your loved one safe and comfortable in your home is your focus. At Med-Mizer, it's our focus, too. When a loved one has an illness or injury, finding the right technology to help keep them at home can make all the difference. We have designed our solutions to provide a more comfortable and safer experience for the user and caregiver.

Please browse through our list of hospital beds and mattresses, and learn how we can provide the right product for your loved one to use in the home. With an industry-leading warranty and products in stock ready to ship, we are confident we can supply you with the right solution for your needs.






NOW YOU DON'T HAVE TO  
CHOOSE BETWEEN A HOME-LIKE  
ENVIRONMENT AND SAFETY,  
QUALITY, AND BUDGET. **BEST OF  
ALL, EVERYTHING COMES FROM  
THE SAME TRUSTED SOURCE  
AND IS MADE IN THE USA.**

Bed Frames . . . . .	pg 4
Bariatric Solutions . . . . .	pg 16
Prevention and Treatment Surfaces . . . . .	pg 18
FlexTilt Chair . . . . .	pg 20
Patient Lift . . . . .	pg 24

With our full line of bed frames we can help your **improve your loved ones** safety and comfort while enhancing **the care you provide.**

**Our Retractable offerings are the choice for space saving design and allow patients to enjoy access to their favorite items while reducing their fall's risk.**

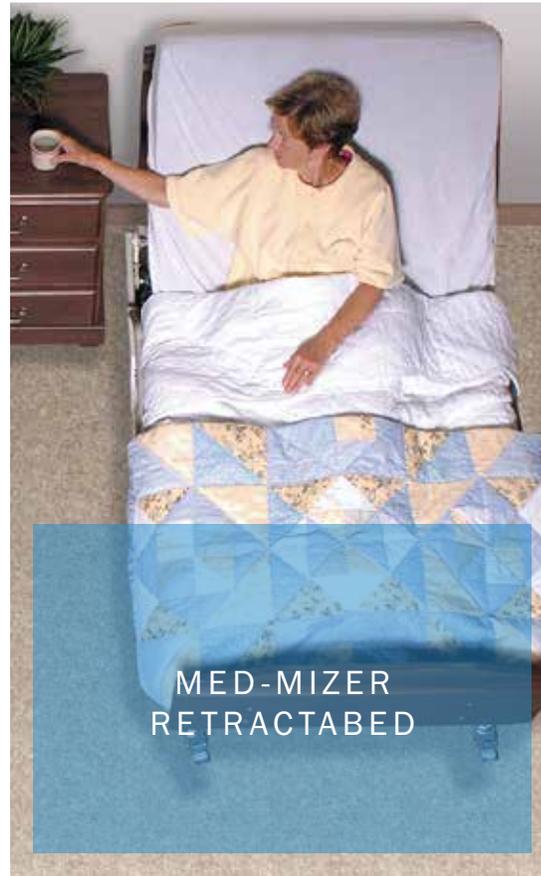

**50%**  
OF FALLS  
OCCUR FALLING  
FROM BED.<sup>(1,2)</sup>

# RETRACTABLE

- Improve safety and satisfaction by placing them at point of living
- Low bed height reduces falls risk
- Clinical Contour® reduces patient migration, alleviates skin shearing, and improves caregiver productivity

## SURELOCK ADDED SAFETY FEATURE AND IMPROVED EASE OF USE FOR LOCKING CASTERS



FIXED




EXPANDABLE

<b>Height</b> .....	9"-25" (w/o surface)
<b>Sleep Surface</b> .....	35/39/42" x 80/84"
<b>Weight Capacity</b> .....	500 lbs (600 lb option)
<b>Options</b> .....	Bed Extender, Battery Back up, SureLock, Trapeze Staff Control, Night Light

1 | Tideisaar, R. The utilization of adjustable low beds in the prevention of falls and injurious falls in long-term care facilities, ECPN, March/April, 2008 30-34.

2 | Oliver D. Bed falls and bed rails-what should we do. Age and Ageing 2002; 31: 415-418 Wide beds: Added inches, added safety <http://www.itlmagazine.com/print/article/wide-beds-added-inches-added-safety> (PDF Download Available). Available from: [https://www.researchgate.net/publication/283151690\\_Wide\\_beds\\_Added\\_inches\\_added\\_safety\\_httpwwwitlmagazinecomprintarticlewide-beds-added-inches-added-safety](https://www.researchgate.net/publication/283151690_Wide_beds_Added_inches_added_safety_httpwwwitlmagazinecomprintarticlewide-beds-added-inches-added-safety) [accessed Mar 14 2018].



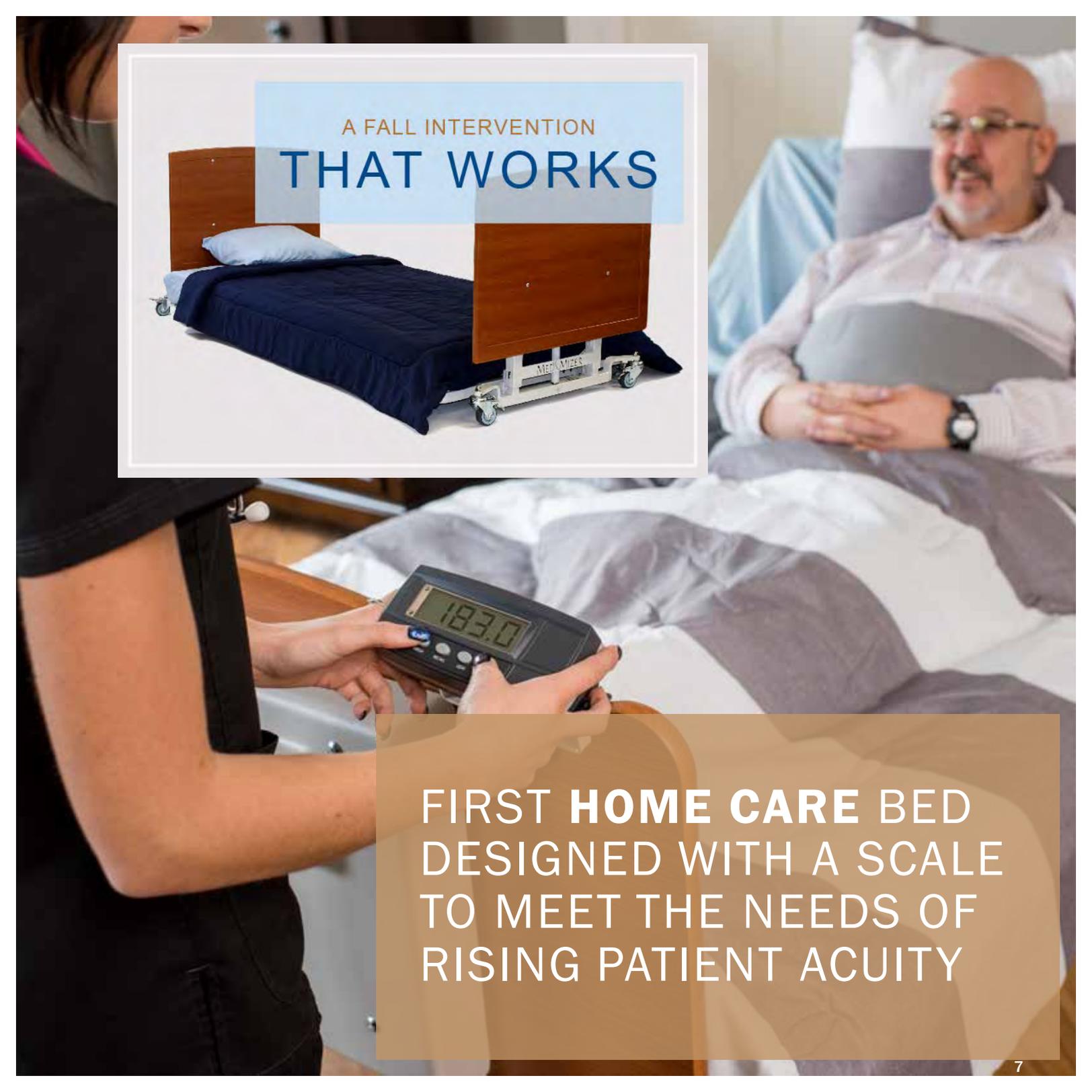

A FALL INTERVENTION  
**THAT WORKS**



**FIRST HOME CARE BED  
DESIGNED WITH A SCALE  
TO MEET THE NEEDS OF  
RISING PATIENT ACUITY**



**Individuals No Longer  
Have to Sacrifice  
Between Looks and  
Functionality**



<b>Height</b> .....	10"-25" (w/o surface)
<b>Sleep Surface</b> .....	35/39/42/48" x 80/84"
<b>Weight Capacity</b> .....	600 lbs
<b>Options</b> .....	Scale, Night Light, Battery Back up



## SELECTCARE BED

- Home-Like Features — enjoy wood aesthetics and hidden casters, and night light to provide your patients the ultimate in "home-like" comfort and amenities
- Ships in 4 pieces for easy assembly into any room in your residence
- Clinical Contour reduces patient migration, alleviates skin shearing, and improves caregiver productivity
- Advanced positioning for added comfort and treatment with One Button Chair, Trend/Reverse Trendelenburg
- Optional scale helps keep accurate patient weight
- Low bed height reduces falls risk
- Comfort Wide expands to 42" without the use of tools to improve patient comfort and dignity



## ActiveCare BED

- **One Button SafeTurn Exit** — exclusive feature that fully automates the bed to turn into a chair and allow the patient to safely enter and exit the bed. Use the built in high/low as a stand assist to offer additional assistance and safety.
- **Full Chair Position** — full chair position to allow the patient the ability to sit in the bed for added comfort. Watch TV, eat a meal, and socialize with family and friends while remaining comfortably in bed.
- **RollBack in Place™** — head of bed can articulate to almost 80 degrees with the assistance of rollback in place to alleviate shearing/friction in the lower back.
- **Home-Like Features** — enjoy wood aesthetics and hidden casters, and night light to provide your patients the ultimate in "home-like" comfort and amenities
- **Advanced "Medical Package" Features** — optional built in scale, out of bed exit alarm/light, ability to sync into nurse call, digital HOB, trendelenburg to create the most advanced "residential looking" hospital bed

### ActiveCare Deluxe

**Height**.....12.75"- 28.75"  
(w/o surface)

**Sleep Surface** .....35" x 80" (AT36)  
42" x 80" (AT42)

**Weight Capacity** .....600 lbs



# Safely Enter and Exit The Bed From a Full Chair Position with the Push of ONE Button





**"This looks like a Rehab professionals' dream bed!...it will be a game changer for independent living within a SNF and along clients using the bed to keep them in their home where they really want to be (comfortable/familiar surroundings)."**

**- RN for SCI Model System**



<b>Height</b> .....	12.5"-27.5" (w/o surface)
<b>Sleep Surface</b> .....	35/42" x 80"
<b>Weight Capacity</b> .....	600 lbs
<b>Options</b> .....	Battery Back up, Trapeze

## ActiveCare Standard Bed


- **One Button SafeTurn Exit** – exclusive feature that fully automates the bed to turn into a chair and allow the patient to safely enter and exit the bed.
- **RollBack in Place™** – head of bed can articulate to 80 degrees with the assistance of rollback in place to alleviate shearing/friction in the lower back.
- **Built in Stand Assist** – use the high/low as a stand assist to set the bed/chair height to an appropriate height to add assistance and safety when getting in and out of bed.
- **Advanced Features** – optional trendelenburg and reverse trendelenburg to help support advanced positioning needs and to help with feet elevation (potentially reduces swelling).

## ActiveCare FixedHeight Bed

- **One Button SafeTurn Exit** — exclusive feature that fully automates the bed to turn into a chair and allow the patient to safely enter and exit the bed.
- **RollBack in Place™** — head of bed can articulate to 80 degrees with the assistance of rollback in place to alleviate shearing/friction in the lower back.

<b>Height</b> .....	12.5"-15.5"-17.5" (w/o surface)
<b>Sleep Surface</b> .....	35/42" x 80"
<b>Weight Capacity</b> .....	600 lbs
<b>Options</b> .....	Battery Back up, 3 and 5" Casters




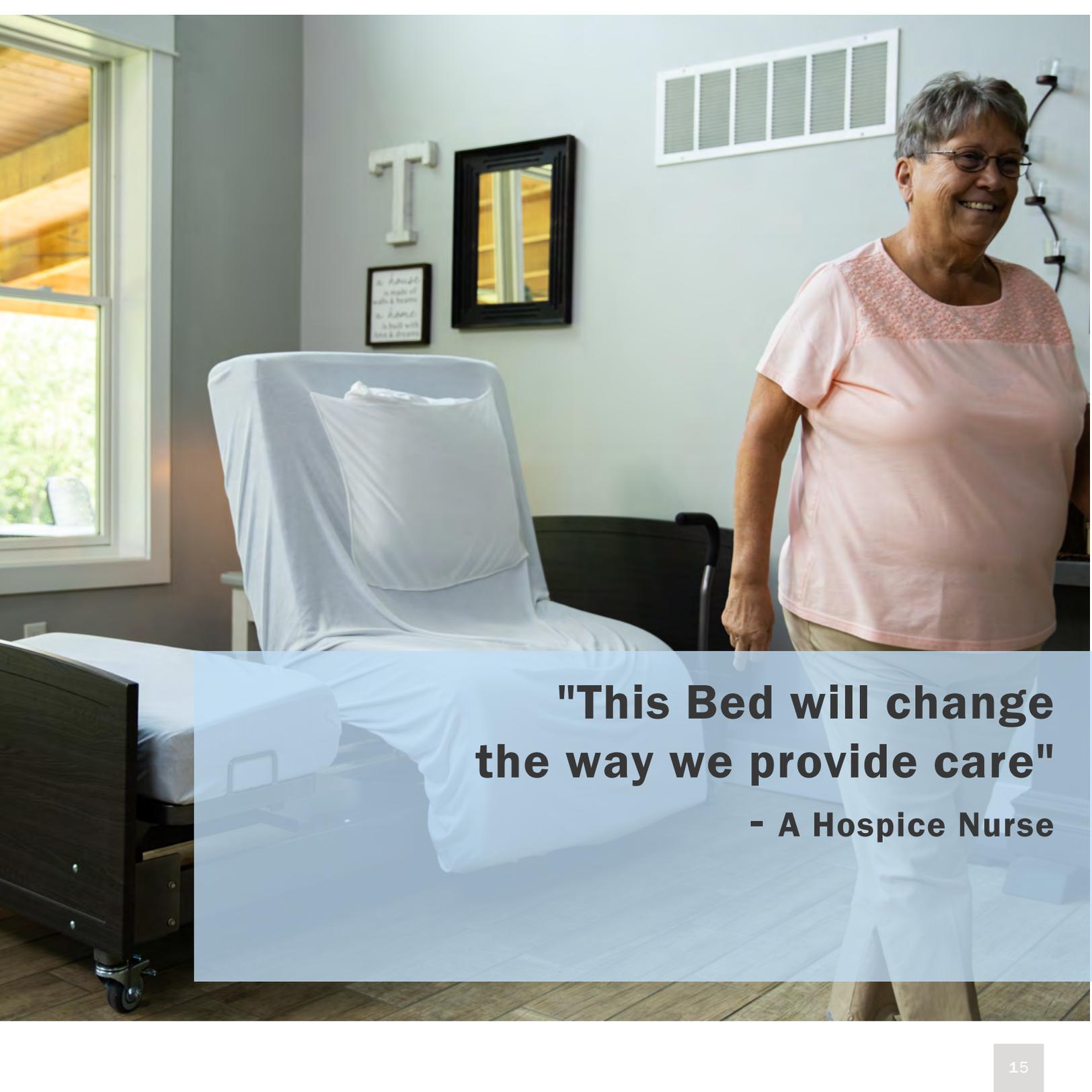

**"This Bed will change  
the way we provide care"  
- A Hospice Nurse**



**OBESITY RATES HAVE  
RISEN 15-25% OVER  
A 10 YEAR PERIOD<sup>(4)</sup>**

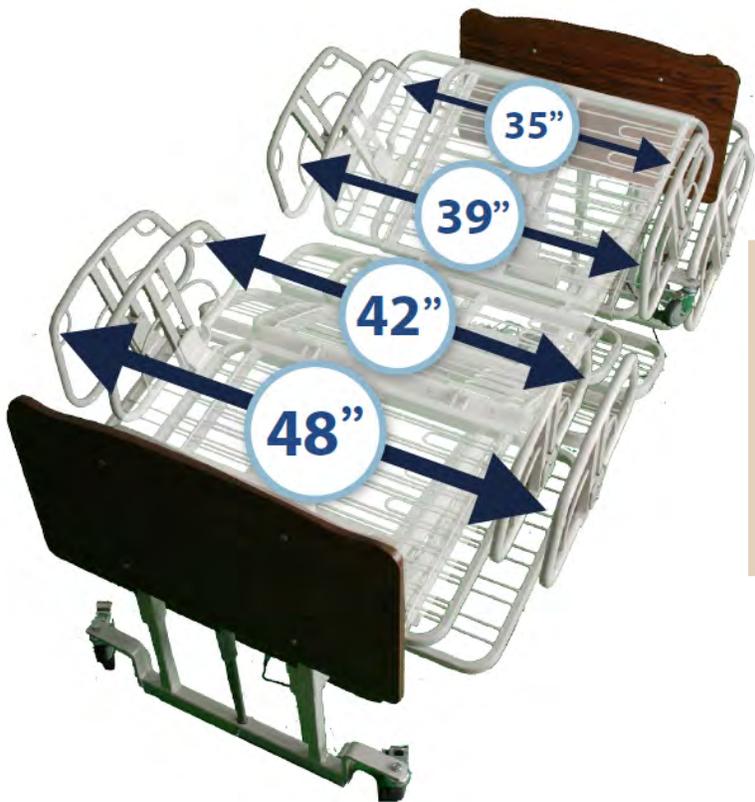


<b>Height</b> .....	9"-25" (w/o surface)
<b>Sleep Surface</b> .....	35/39/42/48" x 80/84"/88"
<b>Weight Capacity</b> .....	800 lbs
<b>Options</b> .....	Scale, Bed Extender, Battery Back up, Trapeze

## COMFORT WIDE EX8000 BARIATRIC BED

- Comfort Wide expands to 48" without the use of tools to improve patient comfort and dignity
- EX8000 ships in 4 pieces for easy assembly into any room in your residence
- Clinical Contour reduces patient migration, alleviates skin shearing, and improves caregiver productivity
- Advanced positioning for added comfort and treatment with Cardiac Chair, Trend/Reverse Trendelenburg
- Optional scale helps keep accurate patient weight
- Low bed height reduces falls risk

# PREVENTION

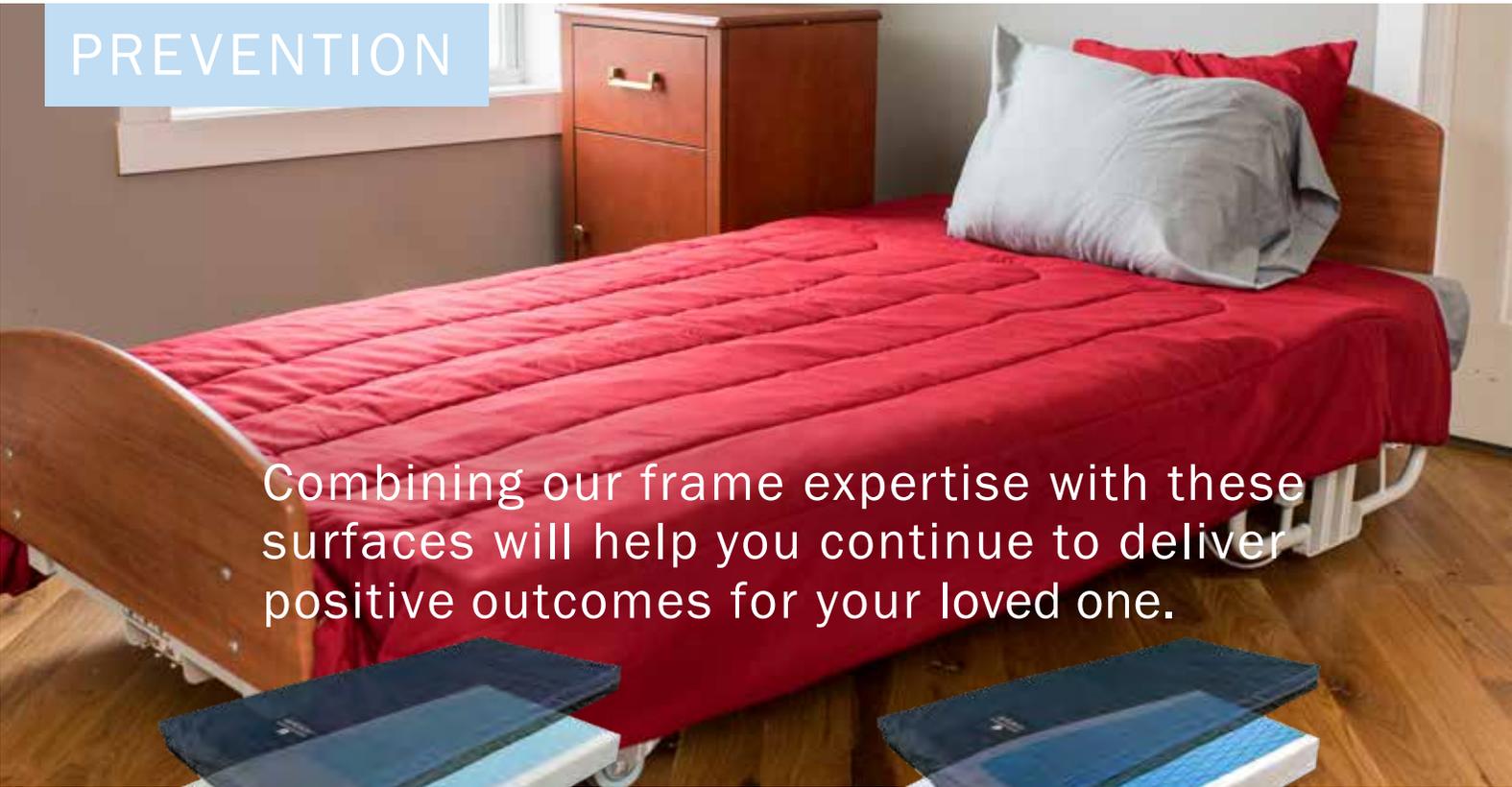

Combining our frame expertise with these surfaces will help you continue to deliver positive outcomes for your loved one.



## Shearless ELITE


- Visco elastic memory foam infused with gel to provide the ultimate in pressure redistribution, enhanced breathability, and to assist with microclimate (heat) management.
- Enhanced heat dissipating material helps minimize perspiration over traditional foam mattresses



## Shearless DELUXE

- Multi-layered pressure redistribution mattress with die-cut top layer featuring high density foam in the torso and head sections.
- Sloped heel section with high density foam offers excellent heel support and pressure redistribution. Additional layer added with channel cut to provide additional pressure relief.

# TREATMENT

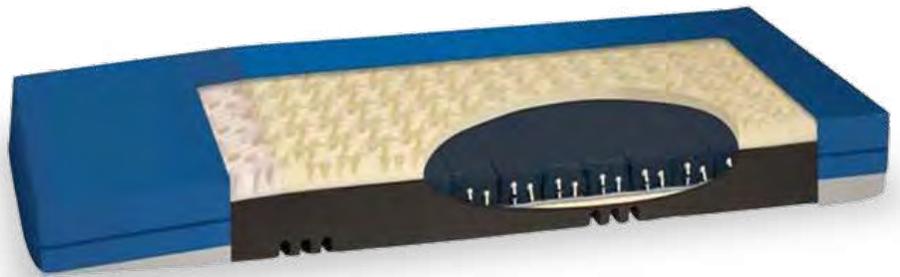

## Shearless StatMatt

- Eight anatomically located horizontal air bladders deliver automatic response to patient movement through open/closed valve controls
- Requires no electricity and is noiseless for added patient comfort
- Firm foam perimeters were designed to enhance transfers and support seating on the edge of the bed.
- Breathable, bi-directional antimicrobial and hypoallergenic stretch coverlet with RF welded seams to prevent fluid ingress.

## Med-Mizer Powered Surface

- Low air loss and alternating pressure with 5 selectable cycle time and CLP
- 8 settings of firmness/softness for comfort
- Auto Firm mode assist with transfers and repositioning
- 20 8" clear "soft" air bladders for added comfort and infection control
- Cell-on-cell technology
- Constant low pressure mode (CLP) for advanced immersion and envelopment
- Auto detection sets appropriate pressure output for patients weight
- Two-way stretch cover with welded seams (anti-fungi, antibacterial, vapor permeable, water resistant)
- 36" (400 lb capacity) and 42" (550 lb capacity)



# WE ARE PROUD TO INTRODUCE THE FLEXTILT CHAIR.

## OUR FLEXTILT CHAIR WILL HELP YOU

- Provide independence and mobility
- Create comfort with unlimited positioning options
- Improve skin/wound care with pressure relief
- Reduce caregiver injuries and efficiency (reduction in transfers, reduced time in repositioning, etc)
- Reduce cost with an affordable chair solution

<b>Frame</b> .....	Gauge steel 11ga - 16ga
<b>Degree of Tilt</b> .....	25°
<b>Degree of Recline</b> .....	30°
<b>Seat Width</b> .....	16", 20", 24"
<b>Seat Depth</b> .....	16.5"-18"
<b>Back Seat Height</b> .....	32"
<b>Seat Height</b> .....	19.5"-21.5"
<b>Weight Capacity</b> .....	350 lbs
<b>Warranty</b> .....	3 year mechanical frame 1 year on components



Ease of use/operation with unlimited positioning for enhanced comfort

# MED-MIZER FLEXTILT CHAIR



**Recline**



**Tilt**



**Tilt + Recline**



"Pressure Relief" strapping

Easy built in handles/control for unlimited positioning

**Vback design for added comfort and support**

Pivoting arms facilitate easier transfers

Recline foot section handle

Built-in length extender for customized seating

Easier maneuverability with steering casters

## MED-MIZER FLEXTILT WC19 TRANSPORT CHAIR



<b>Frame</b> .....	Gauge steel 11ga - 16ga
<b>Degree of Tilt</b> .....	25°
<b>Degree of Recline</b> .....	30°
<b>Seat Width</b> .....	16", 20", 24"
<b>Seat Depth</b> .....	16.5"-18"
<b>Back Seat Height</b> .....	32"
<b>Seat Height</b> .....	19.5"-21.5"
<b>Weight Capacity</b> .....	350 lbs
<b>Warranty</b> .....	3 year mechanical frame 1 year on components

## OUR TRANSPORT CHAIR WILL HELP YOU

- Meets WC19 standards to safely transport patients up to 350 lbs
- Provide independence and mobility
- Create comfort with unlimited positioning options
- Improve skin/wound care with pressure relief
- Reduce caregiver injuries and efficiency (reduction in transfers, reduced time in repositioning, etc)
- Reduce cost with an affordable chair solution





**VBack Design**

Additional Comfort and Support



**Adjustable ArmRest**

Pivoting arms with height adjustability



- Steering Casters – fixed rear casters allow the resident to self propel in a straight direction
- Adjust-ability – customize the chair to the right height and seat depth based on the residents needs
- Independent Leg Rest – optional bolt on accessory allows you to use any standard wheelchair leg rest (no costly proprietary design)
- Made in the SA – an affordable quality product with a commitment to keeping inventory and parts on hand.(no customs or long lead times)

<b>Frame</b> .....	Gauge steel 11ga - 16ga
<b>Degree of Tilt</b> .....	18°
<b>Seat Width</b> .....	16", 20", 24"
<b>Seat Depth</b> .....	16.5"-18"
<b>Back Seat Height</b> .....	32"
<b>Seat Height</b> .....	12.5"-17.5"
<b>Weight Capacity</b> .....	350 lbs
<b>Warranty</b> .....	3 year mechanical frame 1 year on components



## OUR MED-RISER LIFT WILL HELP YOU

- Lift up to 600lbs in a variety of lifting situations (bed/chair transfers, floor, rehab/gait training)
- Create comfort with long lift arm and wider 6-point spreader bar
- Enhanced push handles and push "kick" support bar provides additional leverage
- Powered base and powered height



<b>Safe Working Load .....</b>	<b>600 LBS.</b>
<b>Front Casters.....</b>	<b>4"</b>
<b>Rear Casters.....</b>	<b>4"</b>
<b>Max Height.....</b>	<b>83" (pivot pin)</b>
<b>Low Height.....</b>	<b>19.625" (spreader bar)</b>
<b>Legs Open.....</b>	<b>44" (center of caster)</b>
<b>Legs Closed.....</b>	<b>24" (center of caster)</b>
<b>Overall Length.....</b>	<b>50.2"</b>

# MED-MIZER MED-RISER LIFT



Long Lift Arm improves comfort and allows ultra high lift height

Optional 600lb Scale

6-point spreader with enhanced width for extra comfort and support

Built in controls allow the lift to be used anytime hand controls are in-operable

Powered Height and Base adjustments

Push "Kick" support bar provides additional leverage and safety

## Sling Options



**General Purpose Highback**



**Universal Padded Highback**



**Hygiene**



**Ambulation**

# 3 Easy Steps to a Med-Mizer Home Delivery

Med-Mizer receives requests for our solutions from a variety of sources including Discharge Planners, Home Health Agencies, Physicians, Clinicians, Facilities, Caregivers, Family Members, and Individuals themselves.

## Step 1.

Contact Med-Mizer at 877-867-7365 if you are interested in purchasing products.

## Step 2.

Med-Mizer patient advocates are available to discuss product options and assist in the selection of appropriate products.

## Step 3.

Med-Mizer will schedule product delivery.

- Confirm delivery date and time
- Train end user and caregiver on product usage
- Follow-up with end user and caregiver to ensure product satisfaction

## Med-Mizer Home Care *DIFFERENCE*:

- Manufacture all Med-Mizer Solutions in the USA
- Home care beds are safety stocked for a shorter lead time (usually ships out within 1-2 business days from order)
- Industry leading warranty = more time spent caring for your loved one, not the equipment
- Provide free webinars to highlight product solutions as well as properly train the end user
- Offer white glove delivery to ensure proper set up and installation
- Free technical support to help you troubleshoot any product challenges

PROUDLY MADE IN THE U.S.A.

(877) 867-7365

**MED-MIZER.COM**



is at the forefront of innovating medical technology for post acute care. We are focused on partnering with you to help address healthcare challenges in the following key areas—




FALLS  
PREVENTION



BARIATRIC  
CARE



MOBILITY

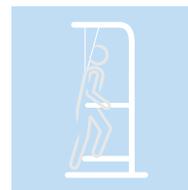

REHABILITATION



WOUND CARE